

GUM SONIC SENSITIVE

SunstarGUM.com

(A)

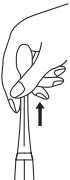(B)

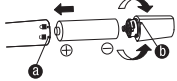(C3)

4101

Made in China

4101LE1020

EN DESIGNED FOR SENSITIVE TEETH AND GUMS. AN ULTRA GENTLE AND EXTRA DEEP CLEAN BETWEEN THE TEETH AND UNDER THE GUMLINE WHERE PLAQUE ACCUMULATES. PROVEN SUPERIOR EFFICACY:

Ultra-thin tip bristles remove plaque 88% deeper in hard to reach areas . Long lasting battery – up to 3 months of use . SONIC battery toothbrush, replaceable head and battery.

A) HOW TO REPLACE THE BRUSH HEAD

Holding base, pull off head. Place new head on base and push down. Replace only with GUM® SONIC refill heads.

B) HOW TO REPLACE THE BATTERY

Turn and pull down handle base. Replace with one new AAA alkaline battery only. When replacing battery, follow symbols on the bottom part for correct positioning. Close the container and dispose the battery according to local regulations.

C) HOW TO CHANGE COLOR RING

1. Remove the brush head. 2. Remove the current silver ring by lifting it up. 3. Select a new ring in your favourite colour and slide it on the handle. Make sure the ring is placed correctly on the handle (see images) and there is no space between parts. 4. Place the head on base and push down.

HOW TO USE

Apply a small amount of toothpaste on the brush head. Insert the brush head inside the mouth. Turn on the brush pressing at the upper central button. Use the same brushing technique as with a manual toothbrush without applying any additional force.

HOW TO CLEAN

Turn off the brush. Rinse brush head with water after brushing.

It is recommended to change the brush head every 3 months.

WARNING. Keep reach out of children. Choking hazard. Small parts. Not for children under 3 years. For children 3 years and older, use under adult supervision only.

PACKAGING CONTENT: 1 toothbrush handle, 1 head, 1 AAA battery, 1 toothbrush cap, 3 color rings

1 vs. a sensitive manual flat trim toothbrush

2 with provided battery for 2 min brushing session twice per day

Distributed by: Sunstar Europe S.A., CH-1163 Etoy
Sunstar Iberia SLU, Calle Aragón 208, 08011 Barcelona, ES

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Jakelija: Sunstar Sverige AB., SE-436 33 Askim
Dystrybuowana przez: Sunstar Europe S.A. Oddział w Polsce, PL- 02-222, Warszawa